

Certificate of Analysis

PEP ME a.s.
PEP ME a.s.
V zářezu 902/4
158 00 Praha-Jinonice
Czechia

Liquilabs s.r.o.
Inovační 122
252 41 Zlatníky-Hodkovice
Czechia
www.liquilabs.cz



[Verify Results Online](#)

Sample Identification

Sample Name Ipamorellin 10 mg
Batch Number
Date Published 2026-08-06 09:09

Results for LYO-0490

Peptides	Result	Unit	Uncertainty	Acceptable Range
Ipamorelin Assay Peptide Screening 0.1% TFA	8.44	mg	[± 0.04]	
Ipamorelin Purity Peptide Screening 0.1% TFA	99.8	%	[± 0.5]	
Ipamorelin Identification by Spectrum (FTIR) Peptide Screening 0.1% TFA	996		[± 5]	
Ipamorelin Identification by RT Peptide Screening 0.1% TFA	0.986		[± 0.005]	
Mass Spectrometry	Result	Unit	Uncertainty	Acceptable Range
Molecular Ion Mass Identification (MS Deconvolution) Mass Spectrometry Identity	711	g/mol	[± 1]	

	Method Specification	
Determination of identity, content and purity of Ipamorelin		
<i>Document number</i> IPA_006_2026	<i>Superseded document</i> -	<i>Number of pages</i> 4

1. Content Assessment

1.1. Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

1.2. Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	280nm


Gradient Program		
Time [min]	A [%]	B [%]
1.5	95	5
13	45	55
13.5	1	99
14.5	1	99
14.51	95	5
16	end	

1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

1.4. Calibration curve

Calibration curve detail	
Quantitative method	External Standard
Calibration Type	Linear
Number of calibration points	5
Force through Zero	Enabled
Weighting Method	None

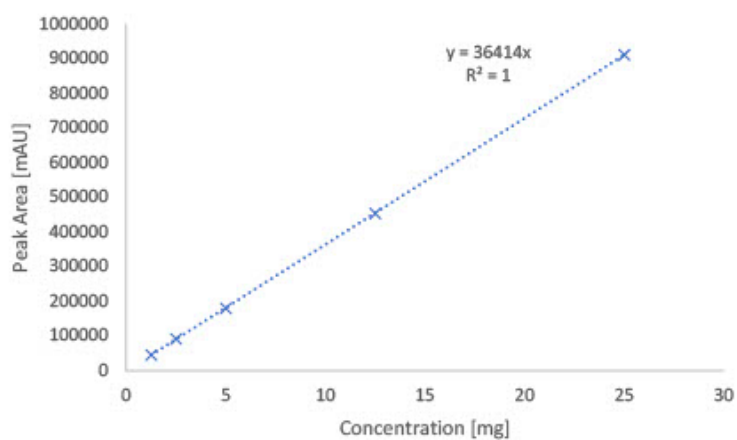

2. Purity assessment

2.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

2.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	225nm

Gradient Program		
Time [min]	A [%]	B [%]
1.5	95	5
13	45	55
13.5	1	99
14.5	1	99
14.51	95	5
16	end	

2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 214 nm.

3. Identity Assessment

3.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

3.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Mass spectrometry	Scan: positive 280-2000 Da

Gradient Program		
Time [min]	A [%]	B [%]
1.5	95	5
13	45	55
13.5	1	99
14.5	1	99
14.51	95	5
16	end	

3.3 Molecular Ion Mass evaluation

Molecular ion mass was determined by deconvolution of multiply charged ESI-MS spectra to calculate the average neutral (zero-charge) molecular mass by equation:

$$M(\text{neutral}) = (z_i((mz_i) - H)) - ME$$

Where:

mz_i - Measured mass of charged particle

z_i - charge

H - proton mass (1.0076 Da)

ME – mass error

Analysis Report





Analysis of quantity and purity of active ingredient by UHPLC with UV detection

Sample Information

Injection Volume : 2
Data File : LYO-0490_031.lcd
Method File : Peptide screening_V7_Group B.lcm
Date Acquired : 8/5/2026 2:16:31 AM

Chromatogram

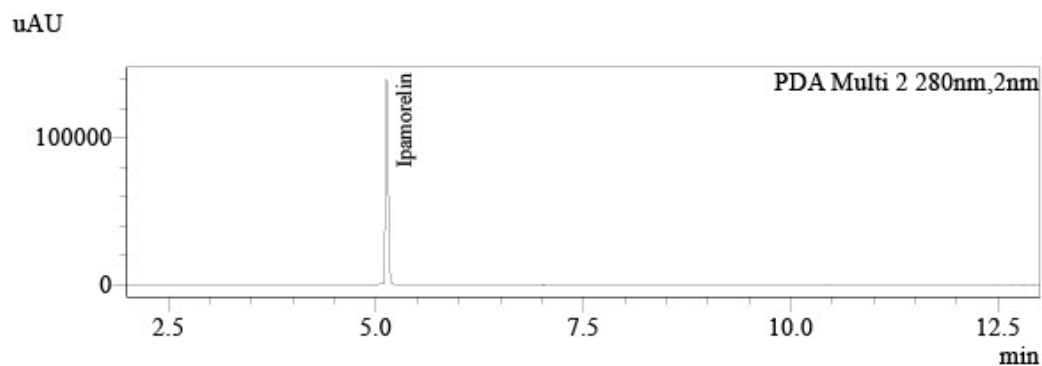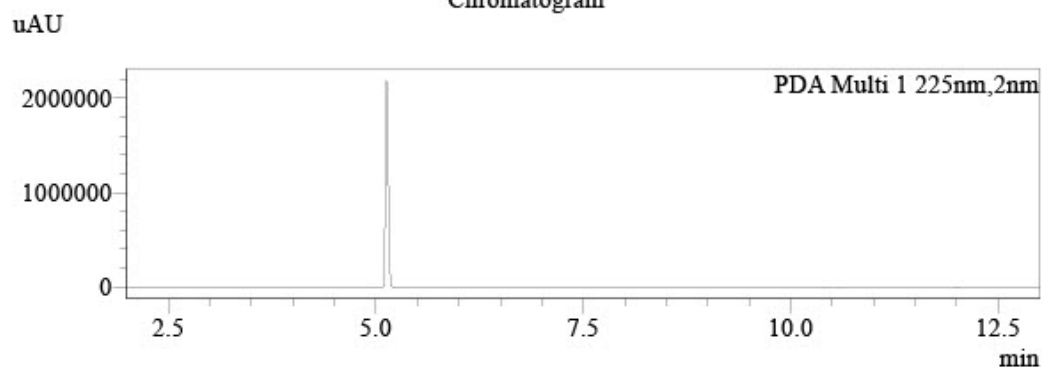

Peak Table

PDA Ch1 225nm					
Name	Ret. Time	Area	Conc.	Unit	Area%
	3.283	1100	0.000		0.021
	4.827	157	0.000		0.003
	4.903	306	0.000		0.006
	4.961	469	0.000		0.009
	5.127	5130521	0.000		99.800
	5.312	1624	0.000		0.032
	5.352	916	0.000		0.018
	5.431	541	0.000		0.011
	5.637	335	0.000		0.007
	5.878	1415	0.000		0.028
	5.912	236	0.000		0.005
	6.163	1496	0.000		0.029
	6.414	743	0.000		0.014
	6.669	419	0.000		0.008
	6.728	349	0.000		0.007
	6.875	164	0.000		0.003
		5140792			100.000

PDA Ch2 280nm

Name	Ret. Time	Area	Conc.	Unit
Ipamorelin	5.127	306972	8.436	mg
		306972		

Analysis Report




LIQUILABS

Analytical services


Analysis of identity of active ingredient by UHPLC with mass spectrometric detection

MS Chromatogram

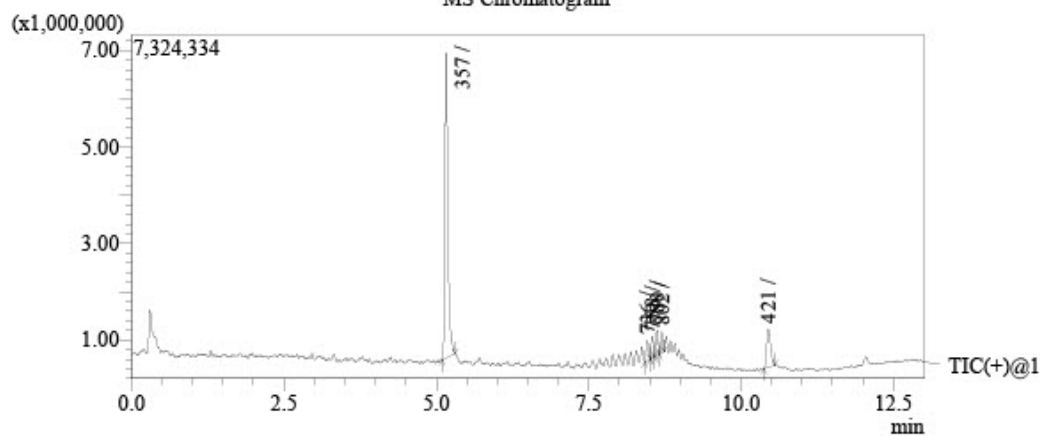


Library Search

Detected compound
Ret. Time: 5.098 Base Peak: 357
Scan Range: 280.00 - 2000.00 Polarity: Positive

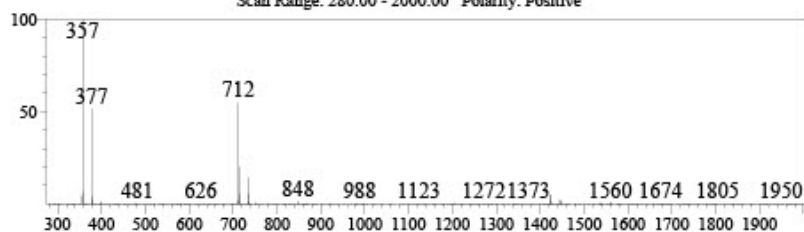

Library: Peptide screening lib
Formula: C₃₈H₄₉N₉O₅ CAS: 170851-70-4 Mol. Weight: 711 Ret. Index: 0
Compound Name: Ipamorelin

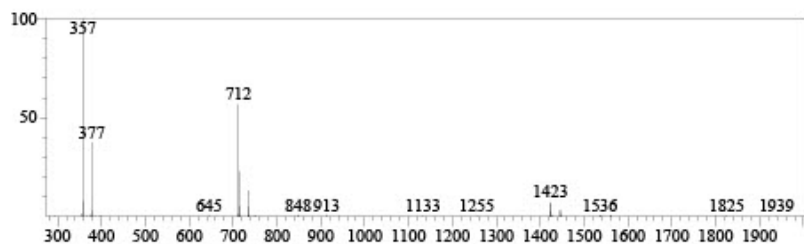

Responsibles



Mr. Ján Galbavý
CEO

Analysis results relate only to the samples tested.

This document shall not be reproduced except in full, without the written approval of Liquilabs s.r.o.