

# Certificate of Analysis

**Nexus Research Laboratories EOOD**

Nexus Research Laboratories EOOD  
Str. Vrabcha No. 17A, Floor 1, Apt. 3,  
Oborishte District  
1000 Sofia  
Bulgaria  
[partnerships@peptralabs.com](mailto:partnerships@peptralabs.com)

**Liquilabs s.r.o.**

Inovační 122  
252 41 Zlatníky-Hodkovice  
Czechia  
[www.liquilabs.cz](http://www.liquilabs.cz)



[Verify Results Online](#)

## Sample Identification

**Sample Name** Tesamorelin 10 mg  
**Batch Number** XD-TSM10-PO05  
**Date Published** 2026-09-19 09:47

## Results for LYO-0887

| Peptides                                                                              | Result | Unit | Uncertainty | Acceptable Range |
|---------------------------------------------------------------------------------------|--------|------|-------------|------------------|
| Tesamorelin Assay<br>Peptide Screening 0.1% TFA                                       | 8.78   | mg   | [± 0.04]    |                  |
| Tesamorelin Purity<br>Peptide Screening 0.1% TFA                                      | 99.8   | %    | [± 0.5]     |                  |
| Tesamorelin Identification by Spectrum<br>Peptide Screening 0.1% TFA                  | 981    |      | [± 5]       |                  |
| Tesamorelin Identification by RT<br>Peptide Screening 0.1% TFA                        | 0.999  |      | [± 0.005]   |                  |
| Mass Spectrometry                                                                     | Result | Unit | Uncertainty | Acceptable Range |
| Molecular Ion Mass Identification<br>(MS Deconvolution)<br>Mass Spectrometry Identity | 5135   | Da   | [± 1]       |                  |

|                                                                                   |                                     |                                 |  |
|-----------------------------------------------------------------------------------|-------------------------------------|---------------------------------|--|
|  | <b>Method Specification</b>         |                                 |  |
| <b>Determination of identity, content and purity of Tesamorelin</b>               |                                     |                                 |  |
| <i>Document number</i><br><br>TESA_008_2026                                       | <i>Superseded document</i><br><br>- | <i>Number of pages</i><br><br>4 |  |

## 1. Content Assessment

### 1.1. Instrumentation

| Module            | Name                 | Serial Number |
|-------------------|----------------------|---------------|
| System Controller | Shimadzu CBM-40 Lite | L221226351398 |
| Degassing Unit    | Shimadzu DGU-403     | NA            |
| Pump              | Shimadzu LC-40B XR   | L22146350580  |
| Autosampler       | Shimadzu SIL-40C XR  | L22216351622  |
| Colum Thermostat  | Shimadzu CTO-40S     | L22236351602  |
| PDA Detector      | Shimadzu SPD-M40     | L22276352808  |
| SQ MS Detector    | Shimadzu LCMS-2050   | O12476200760  |

### 1.2. Chromatographic conditions

| Chromatographic conditions |                                                 |
|----------------------------|-------------------------------------------------|
| Eluent A                   | 0.1% DFA in Water (HPLC, Gradient Grade)        |
| Eluent B                   | 0.1% DFA in Acetonitrile (HPLC, Gradient Grade) |
| Flow rate                  | 0.8 mL/min                                      |
| Program                    | Gradient elution                                |
| Injection volume           | 2 µL                                            |
| Colum Temperature          | 60°C                                            |
| Column                     | Waters XSelect CSH C18, 100x2.1mm 2.5µm         |
| Detection wavelength       | 280nm                                           |


| Gradient Program |       |       |
|------------------|-------|-------|
| Time [min]       | A [%] | B [%] |
| 1                | 98    | 2     |
| 13.99            | 35    | 65    |
| 14               | 5     | 95    |
| 15               | 5     | 95    |
| 15.01            | 98    | 2     |
| 19               | end   |       |

### 1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

### 1.4. Calibration curve

| Calibration curve detail     |                   |
|------------------------------|-------------------|
| Quantitative method          | External Standard |
| Calibration Type             | Linear            |
| Number of calibration points | 5                 |
| Force through Zero           | Enabled           |
| Weighting Method             | None              |

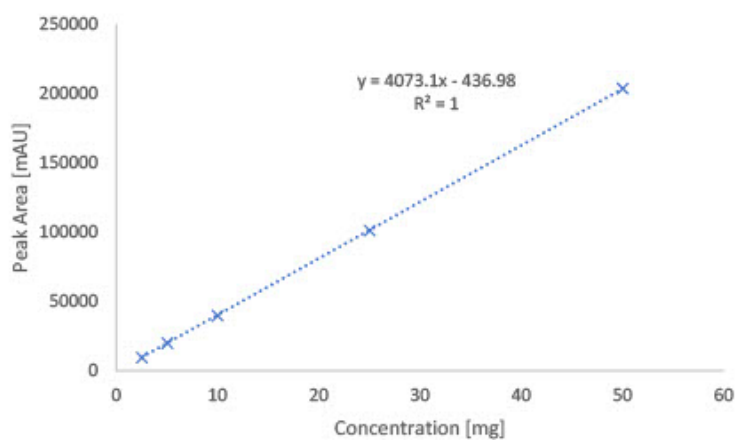

## 2. Purity assessment

### 2.1 Instrumentation

| Module            | Name                 | Serial Number |
|-------------------|----------------------|---------------|
| System Controller | Shimadzu CBM-40 Lite | L221226351398 |
| Degassing Unit    | Shimadzu DGU-403     | NA            |
| Pump              | Shimadzu LC-40B XR   | L22146350580  |
| Autosampler       | Shimadzu SIL-40C XR  | L22216351622  |
| Colum Thermostat  | Shimadzu CTO-40S     | L22236351602  |
| PDA Detector      | Shimadzu SPD-M40     | L22276352808  |
| SQ MS Detector    | Shimadzu LCMS-2050   | O12476200760  |

### 2.2 Chromatographic conditions

| Chromatographic conditions |                                                 |
|----------------------------|-------------------------------------------------|
| Eluent A                   | 0.1% DFA in Water (HPLC, Gradient Grade)        |
| Eluent B                   | 0.1% DFA in Acetonitrile (HPLC, Gradient Grade) |
| Flow rate                  | 0.8 mL/min                                      |
| Program                    | Gradient elution                                |
| Injection volume           | 2 µL                                            |
| Colum Temperature          | 60°C                                            |
| Column                     | Waters XSelect CSH C18, 100x2.1mm 2.5µm         |
| Detection wavelength       | 225nm                                           |

| Gradient Program |       |       |
|------------------|-------|-------|
| Time [min]       | A [%] | B [%] |
| 1                | 98    | 2     |
| 13.99            | 35    | 65    |
| 14               | 5     | 95    |
| 15               | 5     | 95    |
| 15.01            | 98    | 2     |
| 19               | end   |       |

### 2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 225 nm.

### 3. Identity Assessment

#### 3.1 Instrumentation

| Module            | Name                 | Serial Number |
|-------------------|----------------------|---------------|
| System Controller | Shimadzu CBM-40 Lite | L221226351398 |
| Degassing Unit    | Shimadzu DGU-403     | NA            |
| Pump              | Shimadzu LC-40B XR   | L22146350580  |
| Autosampler       | Shimadzu SIL-40C XR  | L22216351622  |
| Column Thermostat | Shimadzu CTO-40S     | L22236351602  |
| PDA Detector      | Shimadzu SPD-M40     | L22276352808  |
| SQ MS Detector    | Shimadzu LCMS-2050   | O12476200760  |

#### 3.2 Chromatographic conditions

| Chromatographic conditions |                                                 |
|----------------------------|-------------------------------------------------|
| Eluent A                   | 0.1% DFA in Water (HPLC, Gradient Grade)        |
| Eluent B                   | 0.1% DFA in Acetonitrile (HPLC, Gradient Grade) |
| Flow rate                  | 0.8 mL/min                                      |
| Program                    | Gradient elution                                |
| Injection volume           | 2 µL                                            |
| Column Temperature         | 60°C                                            |
| Column                     | Waters XSelect CSH C18, 100x2.1mm 2.5µm         |
| Detection wavelength       | 280nm                                           |

| Gradient Program |       |       |
|------------------|-------|-------|
| Time [min]       | A [%] | B [%] |
| 1                | 98    | 2     |
| 13.99            | 35    | 65    |
| 14               | 5     | 95    |
| 15               | 5     | 95    |
| 15.01            | 98    | 2     |
| 19               | end   |       |

#### 3.3 Molecular Ion Mass evaluation

Molecular ion mass was determined by deconvolution of multiply charged ESI-MS spectra to calculate the average neutral (zero-charge) molecular mass by equation:

$$M(\text{neutral}) = (z_i((mz_i) - H)) - ME$$

Where:

$mz_i$  - Measured mass of charged particle

$z_i$  - charge

H - proton mass (1.0076 Da)

ME - mass error

# Analysis Report





## Analysis of quantity of active ingredient by UHPLC with MS detection

### Sample Information

Injection Volume : 2  
Data File : LYO-0887\_\_007.lcd  
Method File : Peptide Screening\_DFA\_Group B.lcm  
Date Acquired : 9/14/2026 7:00:35 PM

### Chromatogram

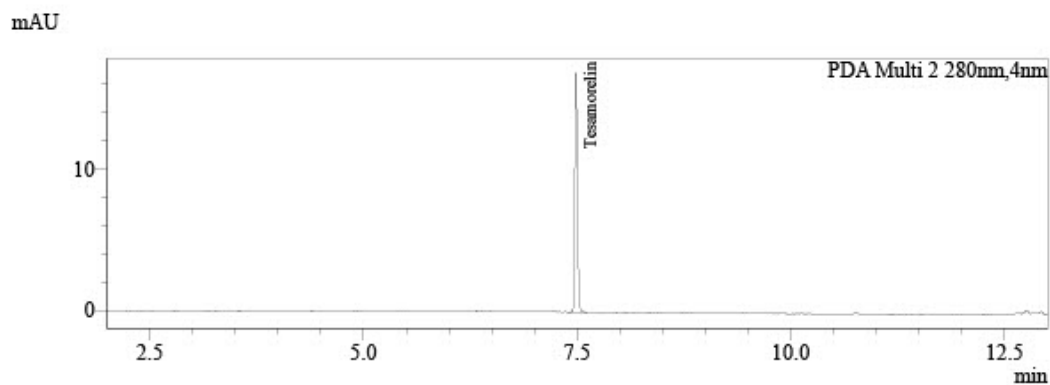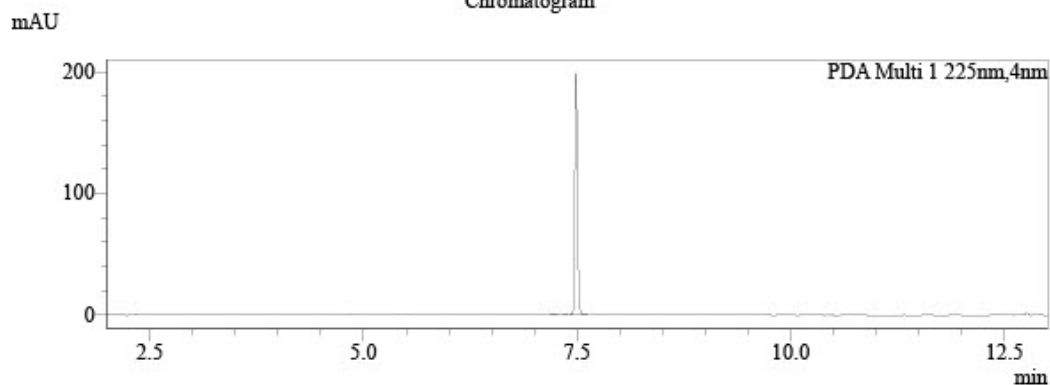

### Peak Table

PDA Ch1 225nm

| Peak# | Name | Ret. Time | Area   | Area%   | Conc. | Unit |
|-------|------|-----------|--------|---------|-------|------|
| 1     |      | 7.220     | 768    | 0.183   | 0.000 |      |
| 2     |      | 7.359     | 165    | 0.039   | 0.000 |      |
| 3     |      | 7.485     | 418457 | 99.778  | 0.000 |      |
| Total |      |           | 419389 | 100.000 |       |      |

### Peak Table


PDA Ch2 280nm

| Peak# | Name        | Ret. Time | Area  | Area%   | Conc. | Unit |
|-------|-------------|-----------|-------|---------|-------|------|
| 1     | Tesamorelin | 7.485     | 35338 | 92.792  | 8.783 | mg   |
| 2     |             | 14.220    | 2745  | 7.208   | 0.000 |      |
| Total |             |           | 38083 | 100.000 |       |      |

# Analysis Report




## Identification of active ingredient by UHPLC Mass Spectrometry

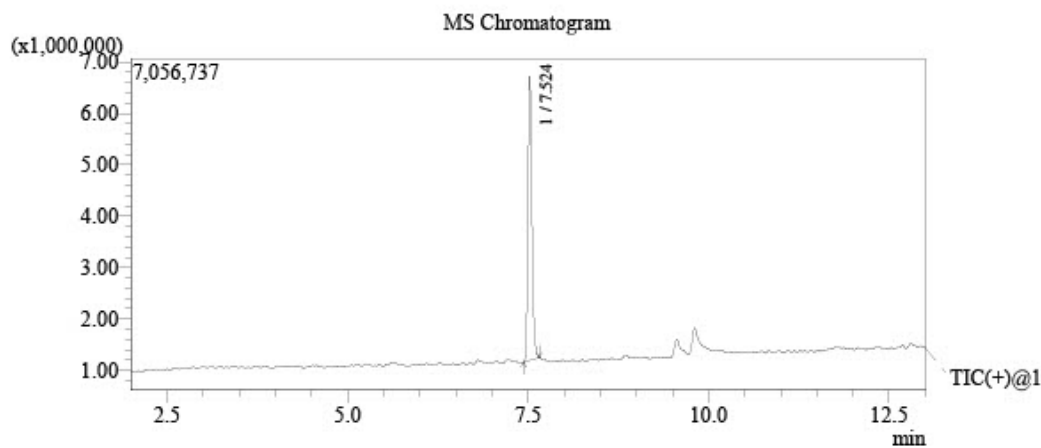


### Library Search

Base Peak: 1047(942743)  
Scan Range: 220.00 - 2000.00 Polarity: Positive

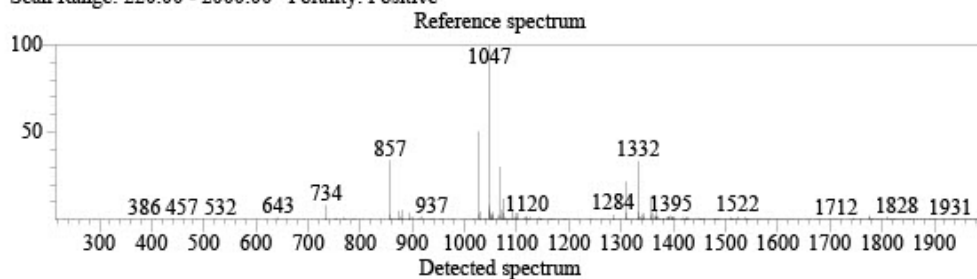

Library: Peptide Screening DFA.lib  
Formula: C<sub>22</sub>H<sub>36</sub>N<sub>7</sub>O<sub>6</sub>S CAS: 218949-48-5 Mol.Weight: 5135(0!=)  
Compound Name: Tesamorelin

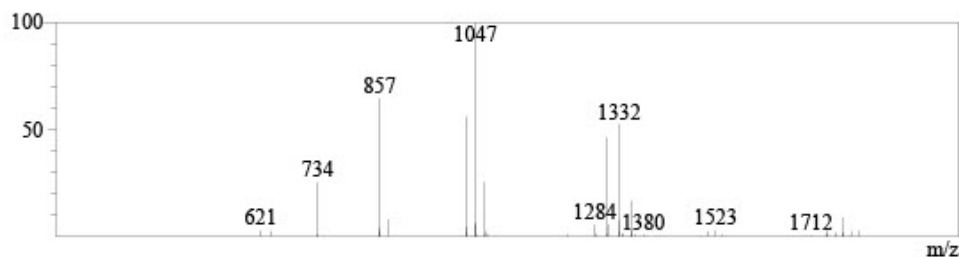

## Responsibles



**Mr. Ján Galbavý**  
*CEO*

Analysis results relate only to the samples tested.

This document shall not be reproduced except in full, without the written approval of Liquilabs s.r.o.