

# Certificate of Analysis

**Peptiva A brand operated by  
Regenevia Pro s.r.o.**

Peptiva A brand operated by Regenevia  
Pro s.r.o.  
Roháčova 145/14  
130 00 Prague  
Czechia  
[support@peptiva.eu](mailto:support@peptiva.eu)

**Liquilabs s.r.o.**

Inovační 122  
252 41 Zlatníky-Hodkovice  
Czechia  
[www.liquilabs.cz](http://www.liquilabs.cz)



[Verify Results Online](#)

## Sample Identification

**Sample Name** MOTS-C 10 mg  
**Batch Number**  
**Date Published** 2026-07-27 09:42

## Results for LYO-0366

| Peptides                                                                                           | Result | Unit  | Uncertainty | Acceptable Range |
|----------------------------------------------------------------------------------------------------|--------|-------|-------------|------------------|
| MOTS-C Assay<br>Peptide Screening 0.1% TFA                                                         | 11.41  | mg    | [± 0.06]    |                  |
| MOTS-C Purity<br>Peptide Screening 0.1% TFA                                                        | 99.5   | %     | [± 0.5]     |                  |
| MOTS-C Identification by<br>Spectrum (FTIR)<br>Peptide Screening 0.1% TFA                          | 986    |       | [± 5]       |                  |
| MOTS-C Identification by RT<br>Peptide Screening 0.1% TFA                                          | 0.998  |       | [± 0.005]   |                  |
| Microbiology                                                                                       | Result | Unit  | Uncertainty | Acceptable Range |
| Bacterial Endotoxin Chromgenic<br>USP<85>/ Eur. Ph. 2.6.14. Bacterial<br>Endotoxin Chromgenic Test | 0.061  | EU/mg | [± 0.001]   | 0 - 0.5          |
| Mass Spectrometry                                                                                  | Result | Unit  | Uncertainty | Acceptable Range |
| Molecular Ion Mass Identification<br>(MS Deconvolution)<br>Mass Spectrometry Identity              | 2174   | Da    | [± 1]       |                  |

|                                                                                   |                                     |                                 |  |
|-----------------------------------------------------------------------------------|-------------------------------------|---------------------------------|--|
|  | <b>Method Specification</b>         |                                 |  |
| <b>Determination of identity, content and purity of MOTS-C</b>                    |                                     |                                 |  |
| <i>Document number</i><br><br>MOTS_006_2026                                       | <i>Superseded document</i><br><br>- | <i>Number of pages</i><br><br>4 |  |

## 1. Content Assessment

### 1.1. Instrumentation

| Module            | Name                 | Serial Number |
|-------------------|----------------------|---------------|
| System Controller | Shimadzu CBM-40 Lite | L221226351398 |
| Degassing Unit    | Shimadzu DGU-403     | NA            |
| Pump              | Shimadzu LC-40B XR   | L22146350580  |
| Autosampler       | Shimadzu SIL-40C XR  | L22216351622  |
| Colum Thermostat  | Shimadzu CTO-40S     | L22236351602  |
| PDA Detector      | Shimadzu SPD-M40     | L22276352808  |
| SQ MS Detector    | Shimadzu LCMS-2050   | O12476200760  |

### 1.2. Chromatographic conditions

| Chromatographic conditions |                                                    |
|----------------------------|----------------------------------------------------|
| Eluent A                   | 0.05% TFA in Water (HPLC, Gradient Grade)          |
| Eluent B                   | 0.0425% TFA in Acetonitrile (HPLC, Gradient Grade) |
| Flow rate                  | 0.9 mL/min                                         |
| Program                    | Gradient elution                                   |
| Injection volume           | 2 µL                                               |
| Colum Temperature          | 55°C                                               |
| Column                     | Waters XSelect CSH C18, 100x2.1mm 2.5µm            |
| Detection wavelength       | 280nm                                              |


| Gradient Program |       |       |
|------------------|-------|-------|
| Time [min]       | A [%] | B [%] |
| 1.5              | 95    | 5     |
| 13               | 45    | 55    |
| 13.5             | 1     | 99    |
| 14.5             | 1     | 99    |
| 14.51            | 95    | 5     |
| 16               | end   |       |

### 1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

### 1.4. Calibration curve

| Calibration curve detail     |                   |
|------------------------------|-------------------|
| Quantitative method          | External Standard |
| Calibration Type             | Linear            |
| Number of calibration points | 5                 |
| Force through Zero           | Enabled           |
| Weighting Method             | None              |

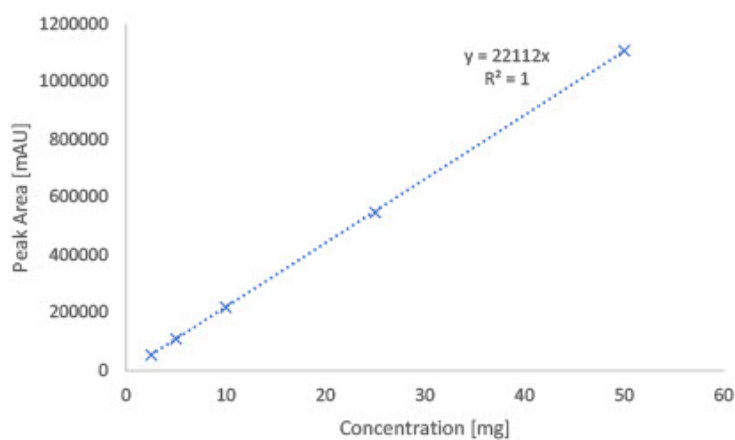

## 2. Purity assessment

### 2.1 Instrumentation

| Module            | Name                 | Serial Number |
|-------------------|----------------------|---------------|
| System Controller | Shimadzu CBM-40 Lite | L221226351398 |
| Degassing Unit    | Shimadzu DGU-403     | NA            |
| Pump              | Shimadzu LC-40B XR   | L22146350580  |
| Autosampler       | Shimadzu SIL-40C XR  | L22216351622  |
| Colum Thermostat  | Shimadzu CTO-40S     | L22236351602  |
| PDA Detector      | Shimadzu SPD-M40     | L22276352808  |
| SQ MS Detector    | Shimadzu LCMS-2050   | O12476200760  |

### 2.2 Chromatographic conditions

| Chromatographic conditions |                                                    |
|----------------------------|----------------------------------------------------|
| Eluent A                   | 0.05% TFA in Water (HPLC, Gradient Grade)          |
| Eluent B                   | 0.0425% TFA in Acetonitrile (HPLC, Gradient Grade) |
| Flow rate                  | 0.9 mL/min                                         |
| Program                    | Gradient elution                                   |
| Injection volume           | 2 µL                                               |
| Colum Temperature          | 55°C                                               |
| Column                     | Waters XSelect CSH C18, 100x2.1mm 2.5µm            |
| Detection wavelength       | 225nm                                              |

| Gradient Program |       |       |
|------------------|-------|-------|
| Time [min]       | A [%] | B [%] |
| 1.5              | 95    | 5     |
| 13               | 45    | 55    |
| 13.5             | 1     | 99    |
| 14.5             | 1     | 99    |
| 14.51            | 95    | 5     |
| 16               | end   |       |

### 2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 214 nm.

### 3. Identity Assessment

#### 3.1 Instrumentation

| Module            | Name                 | Serial Number |
|-------------------|----------------------|---------------|
| System Controller | Shimadzu CBM-40 Lite | L221226351398 |
| Degassing Unit    | Shimadzu DGU-403     | NA            |
| Pump              | Shimadzu LC-40B XR   | L22146350580  |
| Autosampler       | Shimadzu SIL-40C XR  | L22216351622  |
| Column Thermostat | Shimadzu CTO-40S     | L22236351602  |
| PDA Detector      | Shimadzu SPD-M40     | L22276352808  |
| SQ MS Detector    | Shimadzu LCMS-2050   | O12476200760  |

#### 3.2 Chromatographic conditions

| Chromatographic conditions |                                                    |
|----------------------------|----------------------------------------------------|
| Eluent A                   | 0.05% TFA in Water (HPLC, Gradient Grade)          |
| Eluent B                   | 0.0425% TFA in Acetonitrile (HPLC, Gradient Grade) |
| Flow rate                  | 0.9 mL/min                                         |
| Program                    | Gradient elution                                   |
| Injection volume           | 2 µL                                               |
| Column Temperature         | 55°C                                               |
| Column                     | Waters XSelect CSH C18, 100x2.1mm 2.5µm            |
| Mass spectrometry          | Scan: positive 280-2000 Da                         |

| Gradient Program |       |       |
|------------------|-------|-------|
| Time [min]       | A [%] | B [%] |
| 1.5              | 95    | 5     |
| 13               | 45    | 55    |
| 13.5             | 1     | 99    |
| 14.5             | 1     | 99    |
| 14.51            | 95    | 5     |
| 16               | end   |       |

#### 3.3 Molecular Ion Mass evaluation

Molecular ion mass was determined by deconvolution of multiply charged ESI-MS spectra to calculate the average neutral (zero-charge) molecular mass by equation:

$$M(\text{neutral}) = (z_i((mz_i) - H)) - ME$$

Where:

$mz_i$  - Measured mass of charged particle

$z_i$  - charge

H - proton mass (1.0076 Da)

ME – mass error

# Analysis Report



# LIQUILABS

Analytical services


## Analysis of quantity and purity of active ingredient by UHPLC with UV detection

### Sample Information


Injection Volume : 2  
Data File : LYO-0366\_\_006.lcd  
Method File : Peptide screening\_V7\_Group C.lcm  
Date Acquired : 7/20/2026 4:38:06 PM

### Chromatogram

uAU

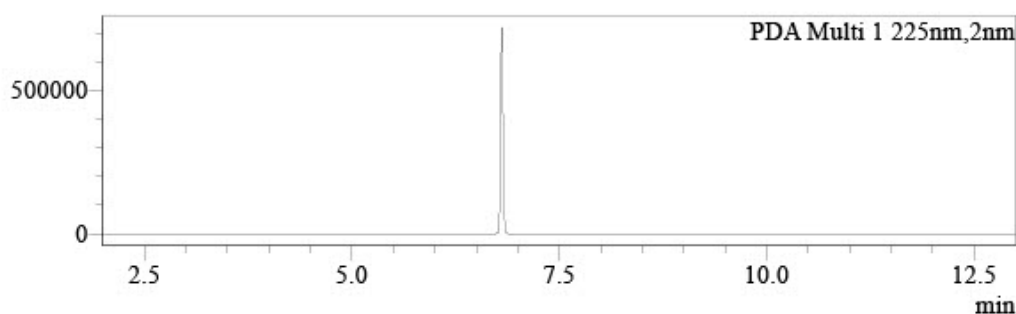

uAU

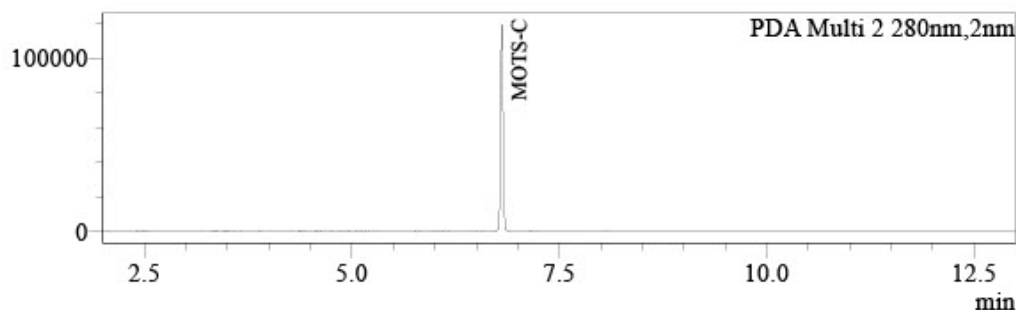

### Peak Table

#### PDA Ch1 225nm

| Name | Ret. Time | Area    | Conc. | Unit | Area%   |
|------|-----------|---------|-------|------|---------|
|      | 5.269     | 929     | 0.000 |      | 0.060   |
|      | 6.185     | 2364    | 0.000 |      | 0.152   |
|      | 6.476     | 1682    | 0.000 |      | 0.108   |
|      | 6.619     | 529     | 0.000 |      | 0.034   |
|      | 6.668     | 160     | 0.000 |      | 0.010   |
|      | 6.806     | 1550954 | 0.000 |      | 99.479  |
|      | 7.217     | 808     | 0.000 |      | 0.052   |
|      | 7.277     | 1174    | 0.000 |      | 0.075   |
|      | 7.620     | 478     | 0.000 |      | 0.031   |
|      |           | 1559078 |       |      | 100.000 |

### Peak Table

#### PDA Ch2 280nm

| Name   | Ret. Time | Area   | Conc.  | Unit |
|--------|-----------|--------|--------|------|
| MOTS-C | 6.806     | 251815 | 11.414 | mg   |
|        |           | 251815 |        |      |

# Analysis Report




**LIQUILABS**

Analytical services


## Analysis of identity of active ingredient by UHPLC with mass spectrometric detection

MS Chromatogram

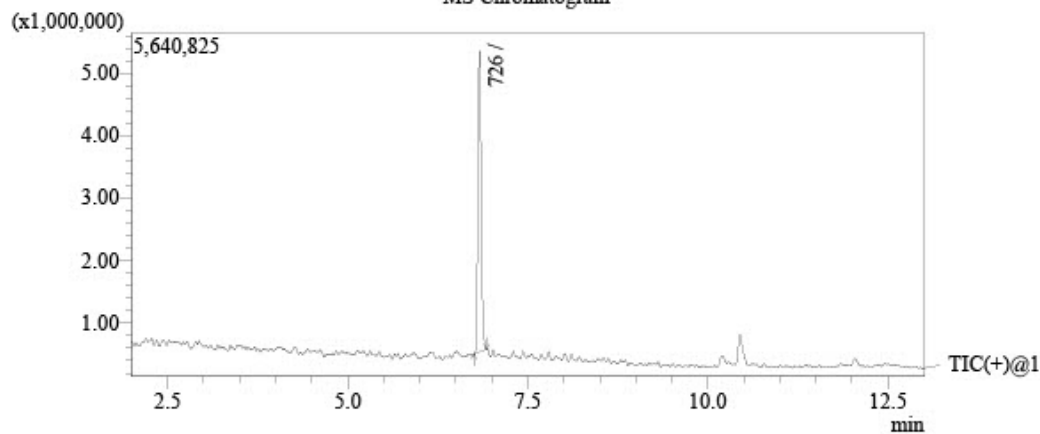


Library Search

Detected compound  
Ret. Time: 6.755 Base Peak: 726  
Scan Range: 280.00 - 2000.00 Polarity: Positive

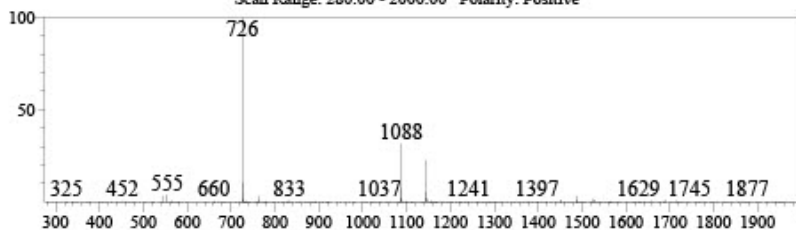

Library: Peptide screening lib  
Formula: C<sub>101</sub>H<sub>152</sub>N<sub>28</sub>O<sub>22</sub>S<sub>2</sub> CAS: 162758-64-6 Mol.Weight: 2174 Ret.Index: 0  
Compound Name: MOTS-C

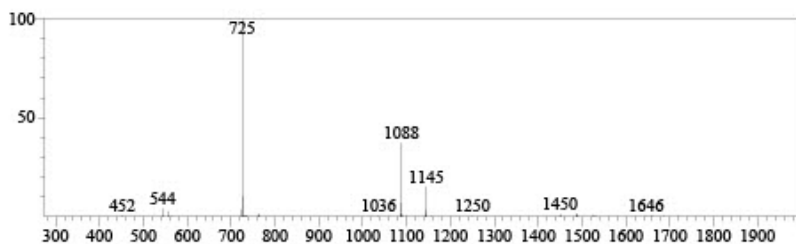

# Endotoxin Determination Report




## File Information

Filename: C:\UVVis-Data\Data\File\_260720\_124259.vqud  
Date/Time: 07/20/2026 12:42:33 PM

## Instrument Information

Instrument Type: UV-1900 Series  
Model (S/N): UV-1900i (A12536253123)

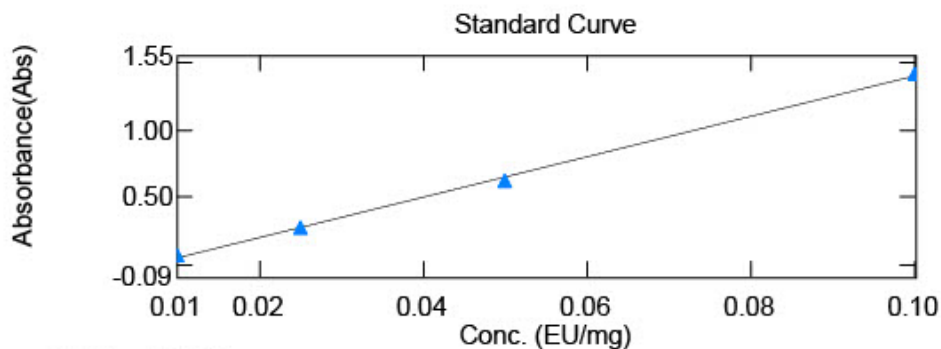

$$y = 15.0474x - 0.104044$$
$$r^2 = 0.99891$$

|    | Sample Name | Conc  | Result | Dilution Rate |
|----|-------------|-------|--------|---------------|
| 1  | LYO-0363    | 0.267 | 0.097  | 20.00         |
| 2  | LYO-0364    | 3.062 | 2.199  | 20.00         |
| 3  | LYO-0365    | 2.603 | 1.855  | 20.00         |
| 4  | LYO-0366    | 0.702 | 0.424  | 20.00         |
| 5  | LYO-0367    | 3.216 | 2.315  | 20.00         |
| 6  | LYO-0368    | 0.487 | 0.263  | 20.00         |
| 7  | LYO-0369    | 0.422 | 0.213  | 20.00         |
| 8  | LYO-0371    | 0.226 | 0.066  | 20.00         |
| 9  | LYO-0372    | 0.243 | 0.079  | 20.00         |
| 10 | LYO-0373    | 0.422 | 0.023  | 50.00         |
| 11 | LYO-0374    | 0.303 | 0.124  | 20.00         |
| 12 | LYO-0375    | 3.200 | 2.304  | 20.00         |
| 13 | LYO-0376    | 0.664 | 0.395  | 20.00         |
| 14 | LIQ-0037    | 0.101 | 0.048  | 10.00         |
| 15 | LYO-0378    | 1.528 | 1.045  | 20.00         |



## Responsibles



**Mr. Ján Galbavý**  
*CEO*

Analysis results relate only to the samples tested.

This document shall not be reproduced except in full, without the written approval of Liquilabs s.r.o.